



GE HealthCare

Revolution Altitude

Rise to a new level of care



Introducing Revolution Altitude

Lives truly change when the advanced becomes ubiquitous.

The demand for highly accessible, efficient CT systems that provide advanced technology to increase workflow, enhance clinical capabilities and provide better patient care has never been greater.

The all-new Revolution™ Altitude is an approachable, customizable introduction to our advanced CT platform. It's been thoughtfully configured for high image quality, increased efficiency and greater clinical capability for those looking to elevate their care day in and day out.



Performance

A thoughtfully designed new imaging chain, combined with deep learning image reconstruction elevates your performance for both advanced and routine imaging.



Efficiency

Smart Workflow solutions provide a higher standard in efficiency through simpler patient set up, positioning, acquisition optimization and post-processing review.



Access

Customizable configurations for those looking to elevate their care day in and day out.

Elevated performance for your practice




Making our revolutionary technology more available means more clinicians can care for more patients faster and with greater ease than ever before.


Which is why Revolution Altitude's imaging chain has been designed to provide greater access to our advanced CT technology featuring the new Lumia detector and our customizable, scalable Performix 2.0 tube.

By combining the new imaging chain with the latest software innovations, including our trusted TrueFidelity DL solution, you will be able to produce high-definition, low-noise CT images with your preferred image texture.

Our new generation Lumia detector highlights image details and delivers consistent image quality across your examinations.

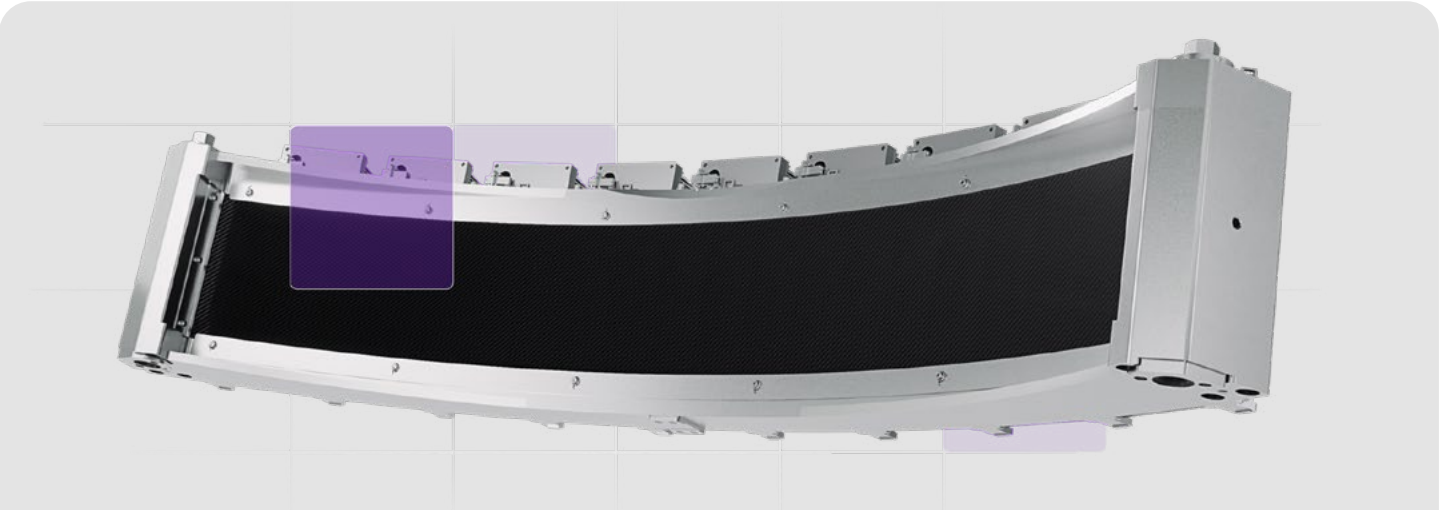

 40 mm of coverage	 0.25 mm spatial resolution ¹	 Up to 30% reduction in electronic noise	 37% reduction in scatter-to-primary ratio for lower artifacts ²
--	--	--	--

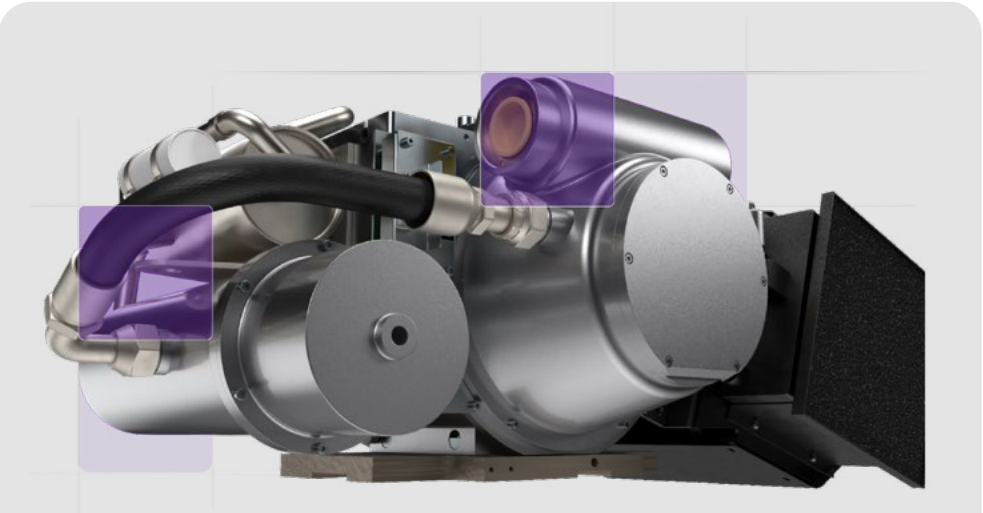

 Up to 0.35 sec rotation speed	 Up to 7 MHU tube capacity
--	--

Our latest Performix 2.0 tube is built for durability and scalability, supporting intensive clinical utilization without compromising performance when scanning a high volume of patients.

Deep learning for a deeper understanding

TrueFidelity DL for true insight

TrueFidelity™ DL is our state-of-the-art AI-based image reconstruction technology that uses a Deep Neural Network (DNN) to generate high-definition, low-noise CT images.

TrueFidelity DL produces images with preferred image sharpness, low-contrast image quality performance and your preferred noise texture, at the same dose.



Filtered Back Projection (FBP)



ASiR-V 40%



TrueFidelity DL-M



Up to **85%** dose reduction at same image quality^{3,4}





Improved low contrast detectability at the same dose⁴



Up to **90%** image noise reduction at the same dose⁵ with preferred noise texture⁷ and boundary sharpness⁸



Significant Contrast Noise Ratio (CNR) improvement at the same dose⁶

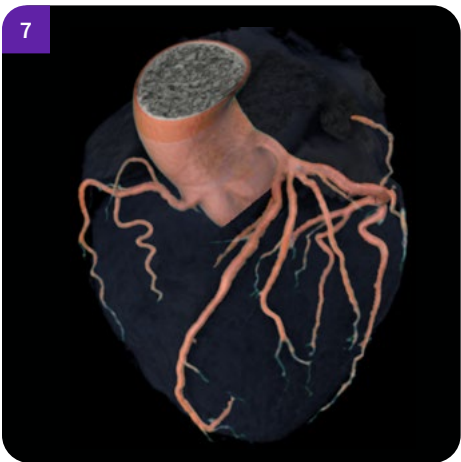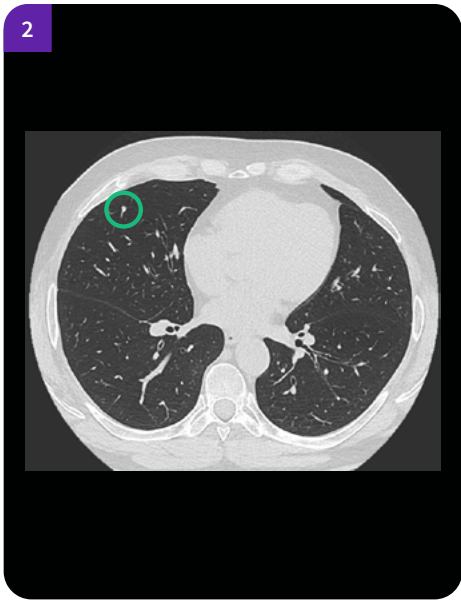

- 1. Bleeding follow up with TrueFidelity DL
- 2. Pulmonary nodule follow up, DLP 81 mGy-cm
- 3. Ruled-out aorta dilatation, coronal view
- 4. Aorta dissection in sagittal view
- 5. 3D volume illumination of lumbar spine stenosis
- 6. Chest, abdominal, pelvis oncological follow up with prostate calcification
- 7. 3D volume illumination of CCTA
- 8. Ankle contusion



A higher standard in efficiency



Efficiency during the scan is critical. The less time it takes to perform the examination, from acquisition to image review, the better it is for the patient and the department. Manual positioning and protocoling can not only add time, it can also lead to sub-optimal image quality and dose.



80% click reduction⁹



56% reduction in positioning time¹⁰



43% increase in table elevation speed¹¹



Up to **308 kg** patient load capacity

Redefine operational efficiency with Smart Workflow

Revolution Altitude's Smart Workflow is redefining patient positioning by combining a gantry-mounted tablet with an AI-powered Xstream Camera.

This integrated approach empowers radiographers with a patient-centric workflow that helps make each scan as comfortable as possible while ensuring the precise dose and consistent imaging you need, all while minimizing system interaction.

Our new Elite table is faster and offers greater flexibility with an increased horizontal range and up to a 308 kg load capacity.



Greater access above all

Revolution Altitude's customizable design allows you to choose the exact premium CT technology that fits your clinical and operational needs today. And you can take comfort in knowing that it will also be ready to grow with your needs tomorrow.



High configurability
to match your
needs and budget
constraints



Choose the technology
you need today,
knowing we'll grow
with you tomorrow







Premium technology
has never been so
accessible





Where detail meets discovery

General imaging

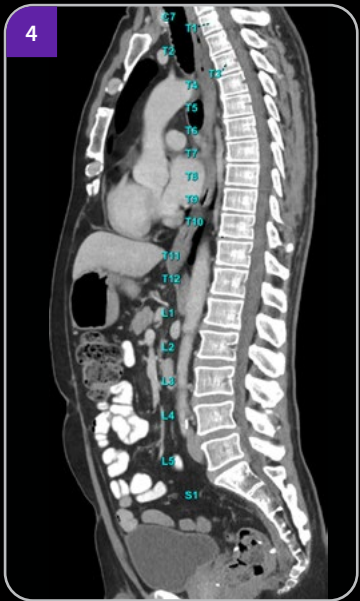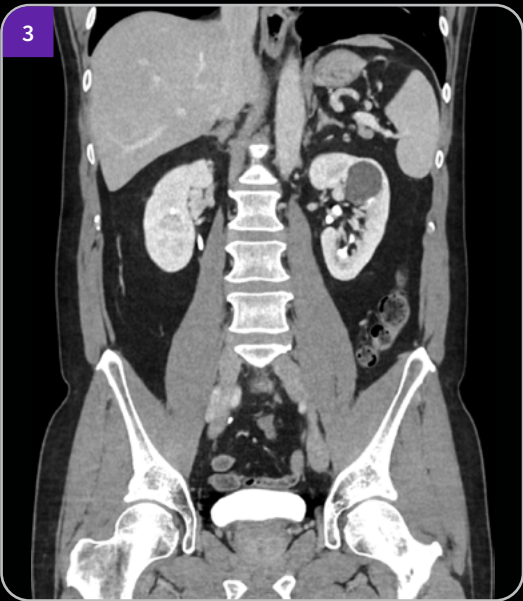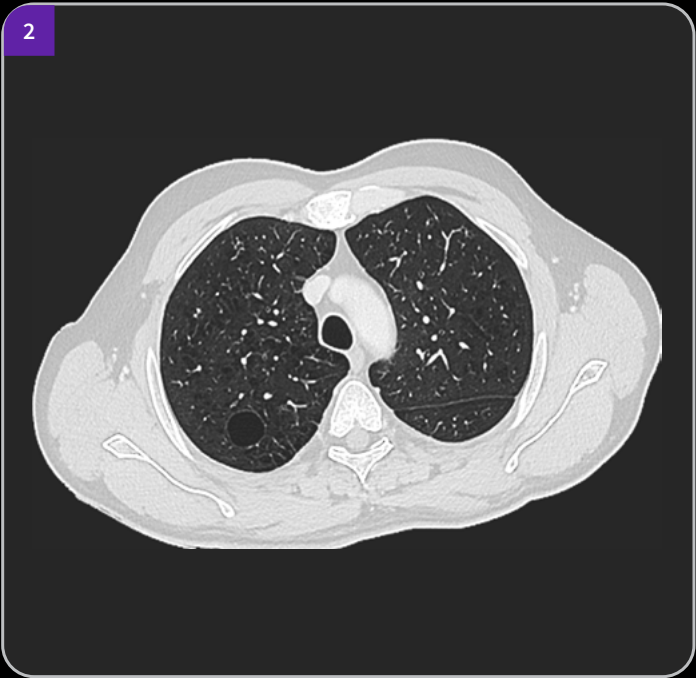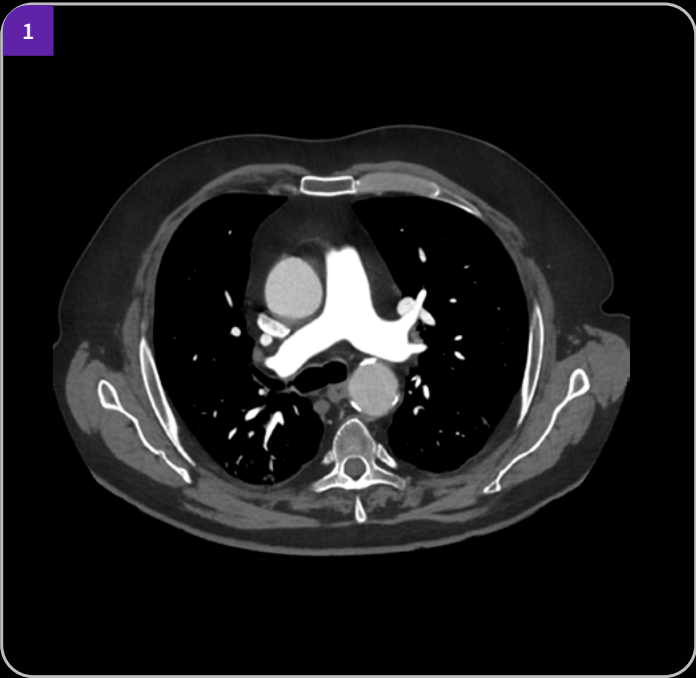

- 1. Thorax exam, 100 kVp
- 2. Oncology follow-up, lung, 100 kVp
- 3. Left renal cyst follow-up
- 4. Chest, abdo, pelvis, onco follow-up




Musculoskeletal

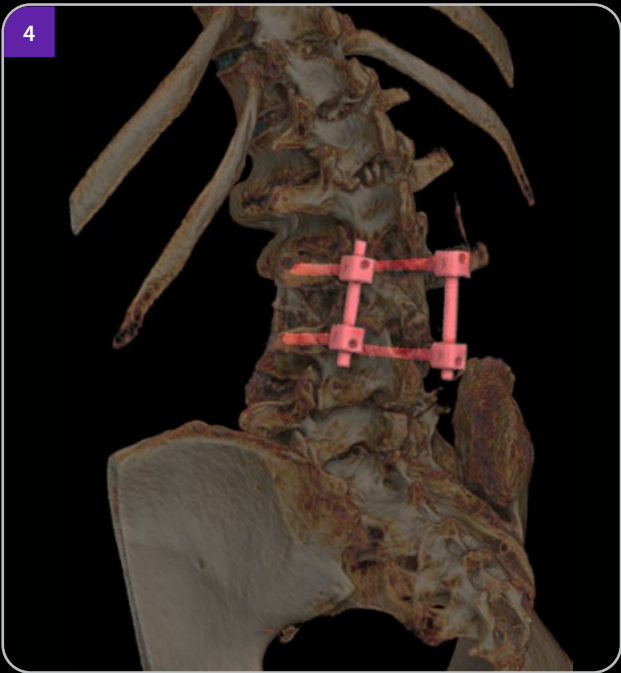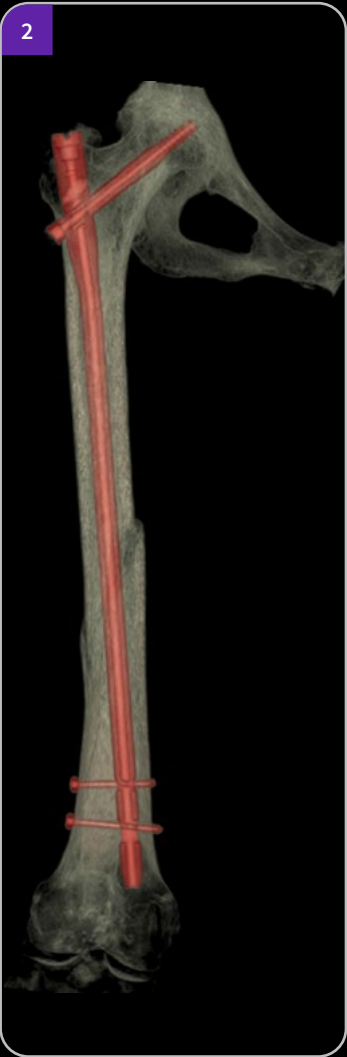

- 1. Sagittal foot, 0.625 mm, 100 kVp
- 2. Femoral bone, 3D VR with MAR
- 3. Lumbar spine stenosis, sagittal view
- 4. Lumbar spine, volume illumination with MAR




Where detail meets discovery

Cardiovascular

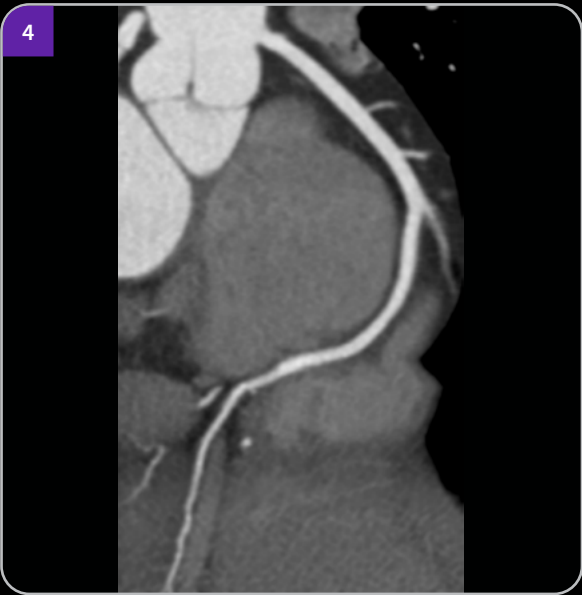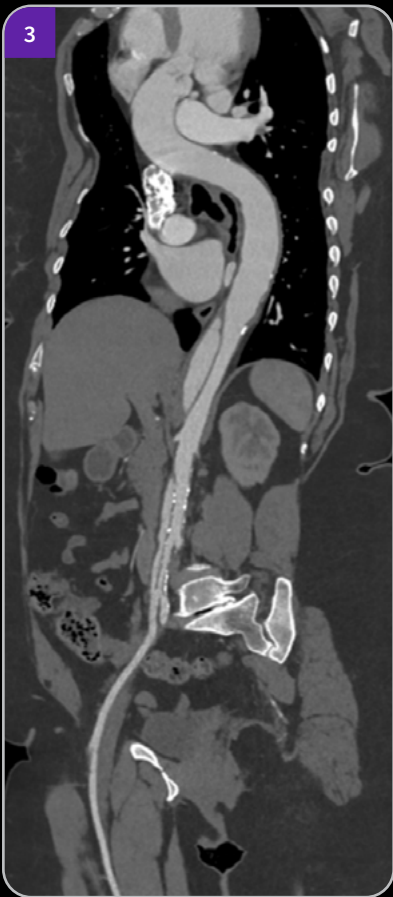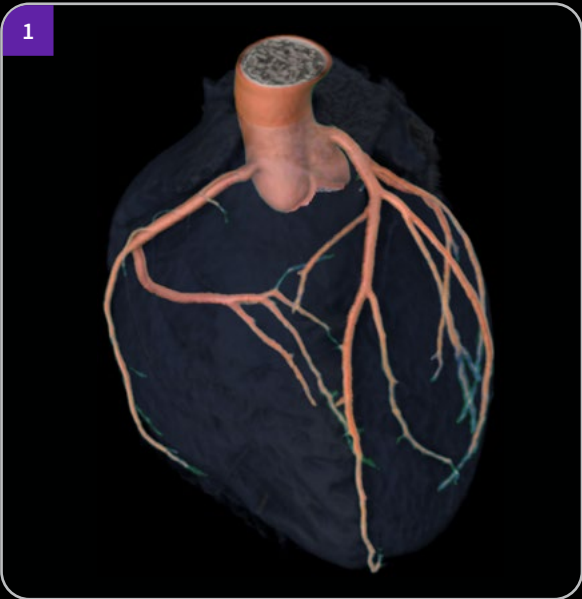

1. Heart, volume illumination, 100 kVp
2. Aorta, volume illumination
3. Aorta, curved view
4. Helical CCTA acquisition, 100 kV

Neuro

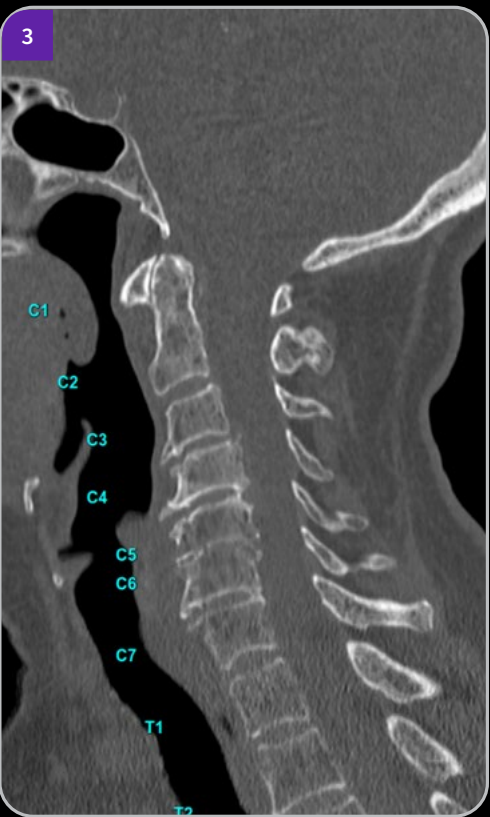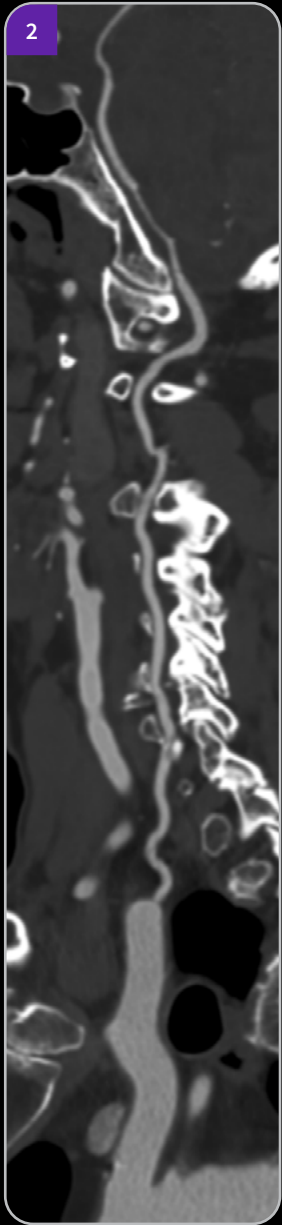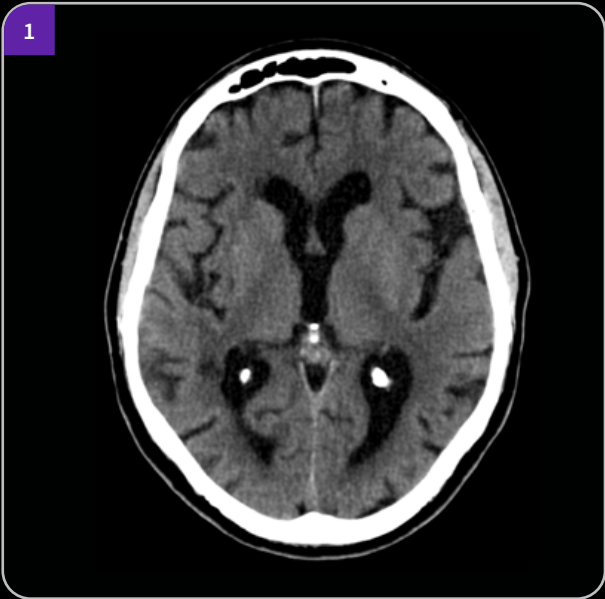

1. Axial brain, 0.625 mm
2. Right vertebral artery, 100 kVp
3. Cervical spine, 0.625 mm
4. Inner ear, 1024 matrix

RISE TO A NEW LEVEL OF CARE



When access to higher-end CT imaging technology is made more obtainable, more lives truly change for the better. This is why Revolution Altitude has been thoughtfully developed with an all-new imaging chain, technology to improve workflow efficiency as well as customizable configurations to fit your budgetary and clinical needs.

Set a higher standard in CT today and rise to a new level of care with Revolution Altitude.

Additional features



Automatic post processing

Auto Views is a suite of automated and AI-based applications, designed to automatically generate anatomically aligned or clinically relevant reformatted views. This enables you to improve image post-processing experience with consistent and reproducible results without manual interactions needed.



All-in-one clinical versatility

Revolution Altitude is built for routine clinical needs yet offers sophisticated clinical capability for cardiovascular, interventional, pediatric and more.



A CT that keeps getting better


CT Smart Subscription ensures your computing platform and software stay up to date, enabling you to keep pace with evolving clinical and workflow innovations through continuous CT upgrades and updates.¹²




Energy saving "green" CT

Energy Saving Mode (ESM) allows the user to place the CT scanner into a low-power mode when not in use. This mode reduces annual energy consumption by up to 47% (in kWh) according to the COCIR SRI 2015 testing procedure with 260 working days per year.

40 mm
detector

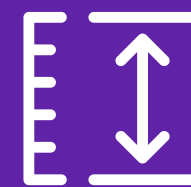

New Imaging
Chain

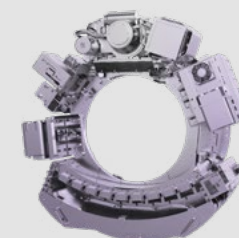

AI-based
Positioning



up to
85%
reduction
in required dose⁴



Up to
.35 sec
rotation speed

Revolution
Altitude


Rise to a new level
of care



Designed for
advanced
procedures



.25 mm
spatial
resolution¹

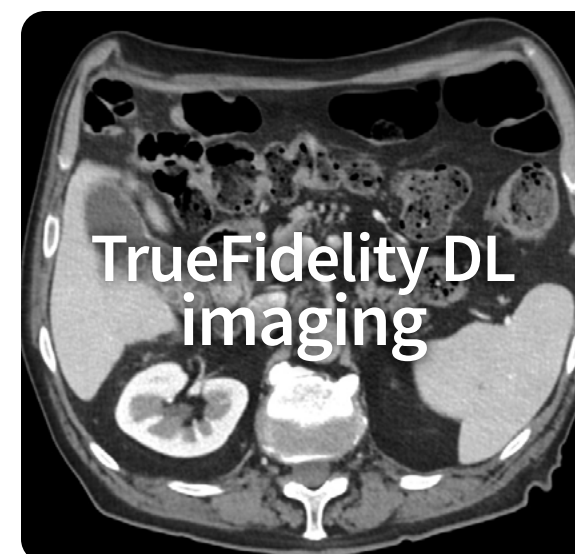

70 cm
Gantry bore with
short geometry

About GE HealthCare Technologies Inc.

GE HealthCare is a trusted partner and leading global healthcare solutions provider, innovating medical technology, pharmaceutical diagnostics, and integrated, cloud-first AI-enabled solutions, services and data analytics. We aim to make hospitals and health systems more efficient, clinicians more effective, therapies more precise, and patients healthier and happier. Serving patients and providers for more than 125 years, GE HealthCare is advancing personalized, connected and compassionate care, while simplifying the patient's journey across care pathways. Together, our Imaging, Advanced Visualization Solutions, Patient Care Solutions and Pharmaceutical Diagnostics businesses help improve patient care from screening and diagnosis to therapy and monitoring. We are a \$19.7 billion business with approximately 53,000 colleagues working to create a world where healthcare has no limits.

Follow us on LinkedIn, X, Facebook, Instagram, and Insights for the latest news, or visit our website <https://www.gehealthcare.com> for more information.

¹ Theoretical calculation results based on MTF 0% value in X/Y.

² Measured with scan conditions of 140 kV, 400 mA, large filter, large focal, 40 mm thickness, and use 48 cm Poly phantom.

³ Image quality as defined by low contrast detectability.

⁴ In clinical practice, the use of DLIR may reduce CT patient dose depending on the clinical task, patient size, anatomical location, and clinical practice. A consultation with a radiologist and a physicist should be made to determine the appropriate dose to obtain diagnostic image quality for the particular clinical task. Low Contrast Detectability (LCD) and Artifact were assessed using reference factory protocols comparing DLIR and FBP. The LCD measured in 0.625 mm slices and tested for both head and body modes using the MITA CT IQ Phantom (CT189, CT191, The Phantom Laboratory), using model observer method.

⁵ In clinical practice, the use of DLIR may reduce CT image noise depending on the clinical task, patient size, anatomical location, and clinical practice. A consultation with a radiologist and a physicist should be made to determine the appropriate dose to obtain diagnostic image quality for the particular clinical task. Image noise was assessed using reference factory protocols comparing DLIR and FBP. The noise reduction was measured in 0.3 mGy CTDI vol and tested for body mode using a 25 X 35 cm OD Oval Body Annulus (part number CTP579, The Phantom Laboratory, Salem, NY) positioned over the uniform section (part number CTP486) of the Catphan® 600 phantom (The Phantom Laboratory, Salem, NY).

⁶ In clinical practice, the use of DLIR may improve CT image contrast-to-noise ratio (CNR) depending on the clinical task, patient size, anatomical location, and clinical practice. A consultation with a radiologist and a physicist should be made to determine the appropriate dose to obtain diagnostic image quality for the particular clinical task. The CNR was assessed using reference factory protocols comparing DLIR and FBP. The contrast-to-noise ratio (CNR) was measured by using the ACR CT accreditation phantom (Gammex 464). Scans were done at 40 mm collimation and 9.2 mGy using an adult abdominal axial scan protocol, which is recommended for clinical use with setting (200 mA tube current, 120 kVp peak voltage, 0.5 s exposure and revolution time, 5 mm reconstruction slice thickness).

⁷ DLIR's noise texture rated as improved over ASiR-V* as demonstrated in clinical evaluations on Revolution CT, consisting of 40 cases and 6 physicians, where each case was reconstructed with both DLIR and ASiR-V and evaluated by 3 physicians. In 92% of the reads, DLIR's noise texture was rated better than ASiR-V's. This rating was based on each individual reader's preference.

⁸ DLIR's image sharpness rated as same as or improved over ASiR-V* as demonstrated in clinical evaluations on Revolution CT consisting of 40 cases and 6 physicians, where each case was reconstructed with both DLIR and ASiR-V and evaluated by 3 physicians. In 100% of the reads, DLIR's image sharpness was rated the same as or better than ASiR-V's. This rating was based on each individual reader's preference.

⁹ From GE HealthCare "AI Auto positioning – Technical White Paper (JB02629XX).

¹⁰ The data was based on comparison between GE's legacy products (16 ch and 64 ch scanner) and Revolution Ascend in the three institutions using a pilot product and selected routine head and body. The data set of this comparison was 838 exams for legacy products and 1387 exams for Revolution Ascend. The time saving value may not be effective for all institutions depending on the clinical practice. Defined the scan setting time is from "Open new patient" to "Confirm setting for Scout".

¹¹ Measured for full range motion and compared with previous generation table GT1700.

¹² Software available to customer is dependent on the software package purchased by customer. Availability may differ from country to country.

