

Shanghai United Imaging Healthcare Co., Ltd.
Copyright © Shanghai United Imaging Healthcare Co., Ltd. All Rights Reserved.

Shanghai, China
2258 Chengbei Rd., Jiading District, Shanghai, 201807.

Email | info.global@united-imaging.com
Business Consultation | +86 (21) - 67076666
After-sales Service | 4006 - 866 - 088

Edition ID | 88000057 - MPD - BRE - 01



uCT 960+

**Attainable Intelligence
Simply Masterful**

ABOUT UIH

Shanghai United Imaging Healthcare Co., Ltd. develops and produces a full portfolio of advanced medical imaging and radiotherapy equipment and offers medical IT and intelligent solutions. Founded in 2011 and headquartered in Shanghai, the company has subsidiaries and R&D centers across China, US, and other parts of the world. Half of its 3,300 employees work on research and development. To date, more than 3,000 hospitals worldwide are using UIH products.

To learn more, visit <https://www.united-imaging.com>



uCT 960+

Leading technology platform, the core performance of comprehensive ascension.

16_{cm}

Coverage Z-detector

***0.25_s**

Rotation Speed

uAI

Vision

320

Rows

640

Slices


82_{cm}

Gantry Aperture

0.25s is options


16cm coverage Z-detector

With a 16cm detector high-quality whole-brain, whole-coronary axial scan images can be obtained within a single-rotation scan.

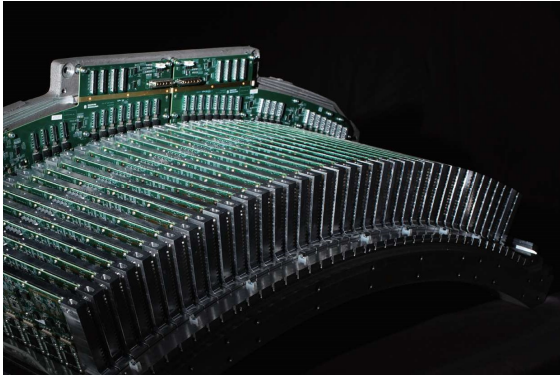

One-beat cardiac imaging and a whole organ perfusion acquisition like a perfusion acquisition of the brain, liver, kidneys, pancreas and other organs.

Volumetric imaging, effective removal of spiral artifacts, enable low-dose imaging with excellent image quality.

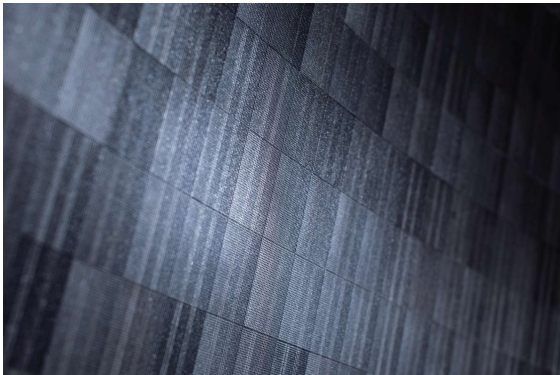

The Real 3D Full reconstruction technology can reduce the cone beam artifact and ensure CT number uniformity across the entire volume of detector coverage.


The 3D ASG (Anti-Scatter Grid) effectively shields the scattered photons without compromising the utilization ratio of non-scattered X-ray.

*0.25s rotation speed

0.25s industry leading rotation speed for high native temporal resolution.



Integrated machine gantry combined with direct drive motor system, showing high performance and high stability.

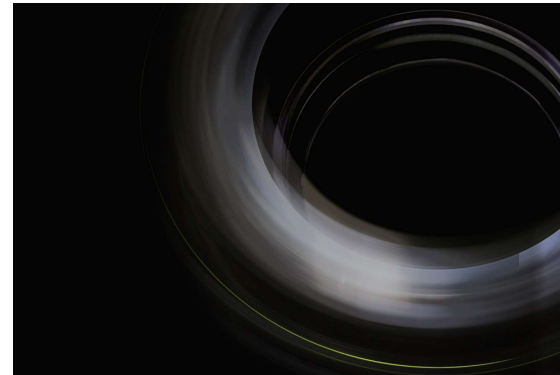

ASIC structure design enables the detector to have fast response and data reading ability, with high sampling rate and processing efficiency.

0.25s is options



uAI Vision:AI-empowered Scan Navigation System

The uCT 960+ with the uAI Vision scan navigation system provides an efficient,standardized and personalized scanning experience for each patient.

The uAI Vision 3D Camera builds real-time digital models for every patient utilizing an AI-based deep learning technology.

Anatomical structures of the patient can be identified with different patient position.

All Clinical Solutions

The uCT 960+ is United Imaging's ultra-premium computed tomography scanner, enabling you to reach a new height in CT imaging. With industry-leading specifications and AI -empowered technologies, the uCT 960+ provides new CT capabilities with impressive image quality, low dose and intelligent applications.

The uCT 960+ provides robust clinical solutions for you and your patients, even during the most challenging clinical scenario period, and seeks successful examination for all patients.

Cardiovascular

Unrestricted cardiac imaging

Stroke

One contrast injection, from anatomical to functional information

Trauma

Fast speed emergency solution

Oncology

Dual energy / Perfusion / 4D dynamic CTA

First-class Experience

Efficient / Minimize the radiation dose / Comfortable



Cardiovascular

The uCT 960+'s gantry supports the industry leading rotation speed of 0.25s. Along with the 16cm detector coverage, the uCT 960+ can image the whole heart in one heart beat.



CardioCapture

Intelligent motion correction with AI-empowered coronary artery extraction for precise movement tracking.



25 ms effective whole-heart temporal resolution

160 mm whole heart coverage with *0.25 s rotation speed to achieve robust, low dose and high quality cardiac imaging within one heart beat.



Intelligent workflow

The uCT 960+ has a streamlined and intelligent workflow, from scanning preparation to image reconstruction, with fully integrated features on the console, including CardioAssist, CardioAdapt, CardioXphase and CardioCapture.



Cardio-vascular imaging

The uCT 960+ allows a mixed acquisition mode of ECG-gated axial scans and non-ECG-gated axial or helical imaging with a fast switching time of less than 2 seconds.


0.25s is options




0.25s is options

One-beat Cardiac Scan

Severe Arrhythmia 52-238bpm
Scan type: Axial one-beat

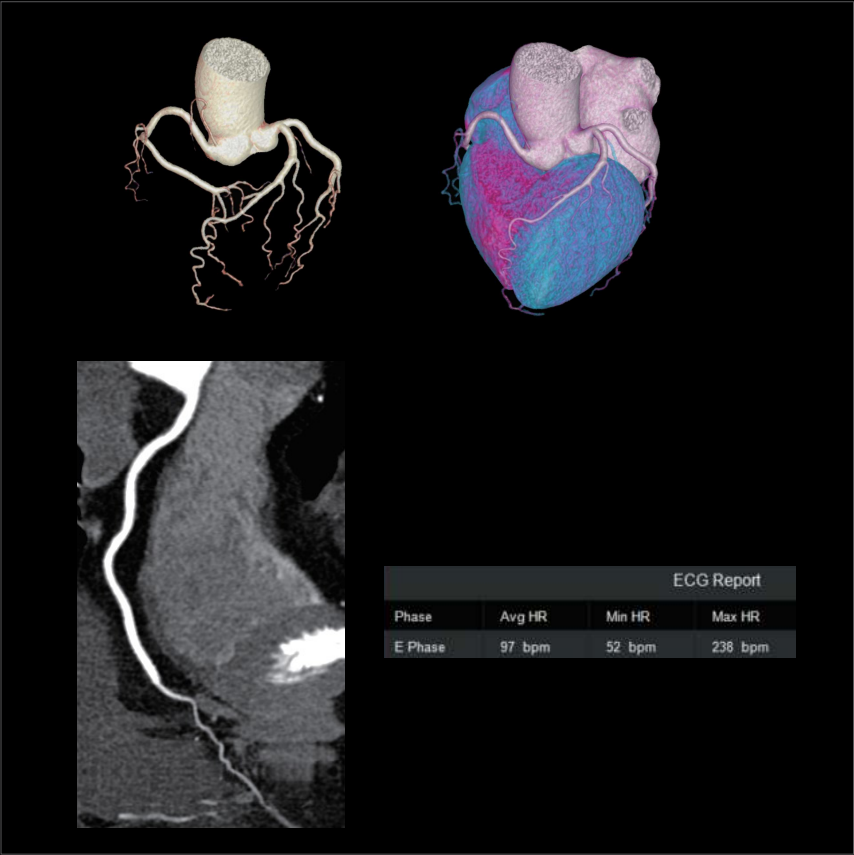

High heart rate with myocardial bridges 99-104bpm
Scan type: Axial one-beat

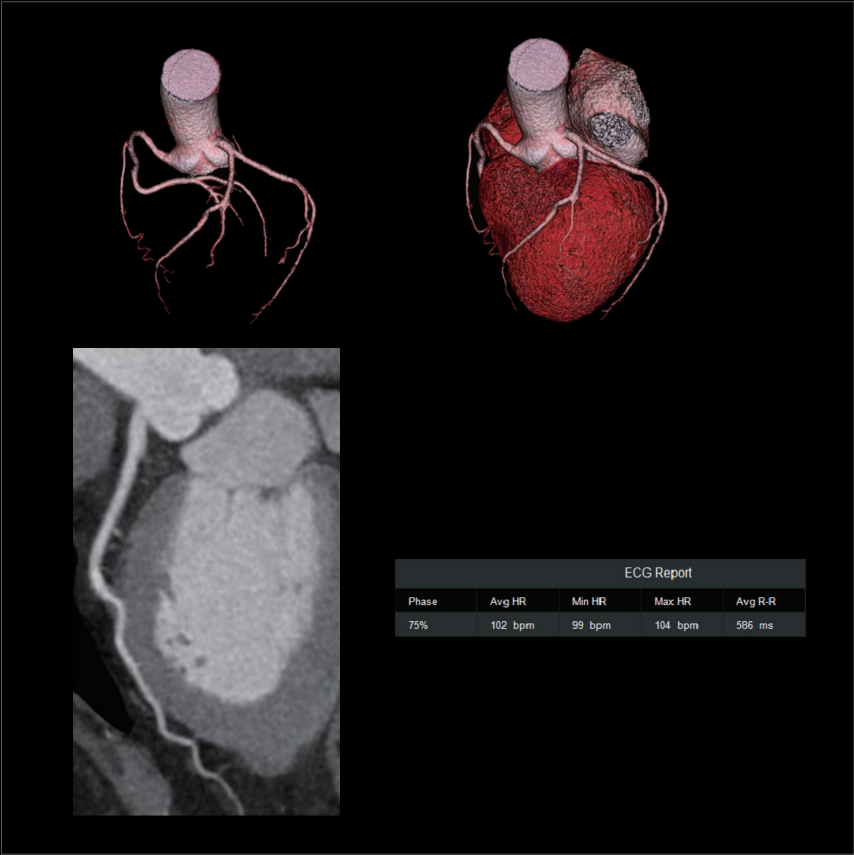

Multi-site Vascular Imaging

One-stop cardio-cerebral vascular Imaging


Multi-site vascular imaging with the use of just a single injection of contrast media. One-stop cardiac-cerebral vascular intelligent analysis post-processing to help physicians diagnose disease with more confidence with a simplified workflow.



Triple Rule Out

One stop acquisition for the evaluation of coronary arteries, pulmonary artery and aorta for triple-rule-out examination.



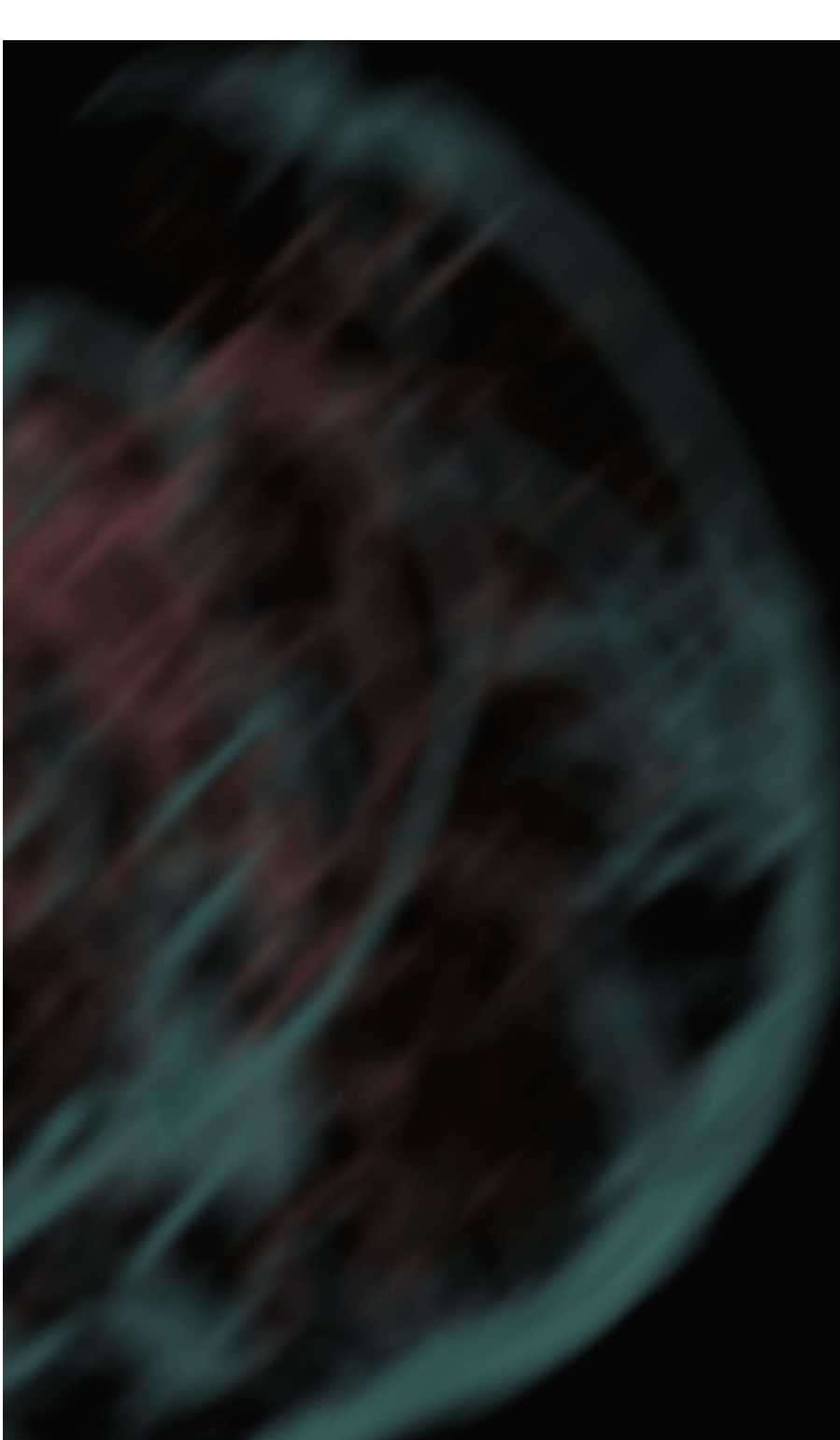

Stroke

One Contrast Injection, From anatomical to functional information
One-stop stroke imaging solutions for physicians to assist with delivering the faster care, sparing as much tissue as possible for the patient.


NCCT

with a 16 cm z-coverage, single scan coverage of the brain. The images can support the evaluation of intracerebral hemorrhage rule-out.




CTA

The CTA images are extracted from the arterial phase of the CTP scan. single scan. single injection.

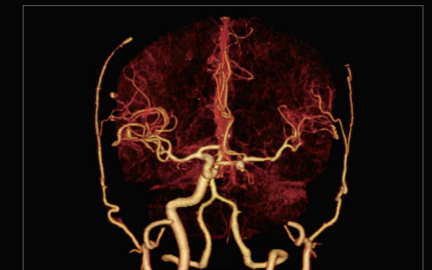


4D Dynamic CTA

Acquires 4D dynamic angiographic image in one injection, auto extraction and separation of artery and vein.

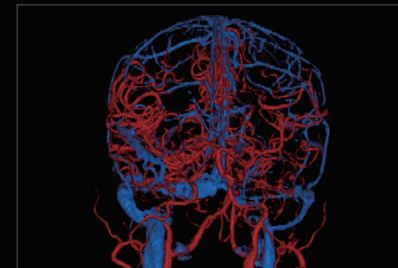

CTP


Automatically calculate CBV, CBF, TTP, MTT, Tmax to help physicians diagnose disease with more confidence.

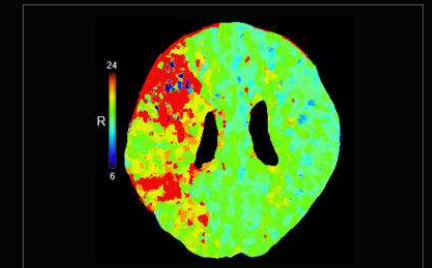

4D one-stop stroke imaging solution


A 4D dynamic vascular and cerebral perfusion examination was performed with one injection

Automatic extraction of arteries and veins, and provide CBV, CBF, MTT, TTP, Tmax perfusion data, support symmetric ROI time density curve and one-click generation of statistical parameters,support CSV format to export ROI statistical parameters.

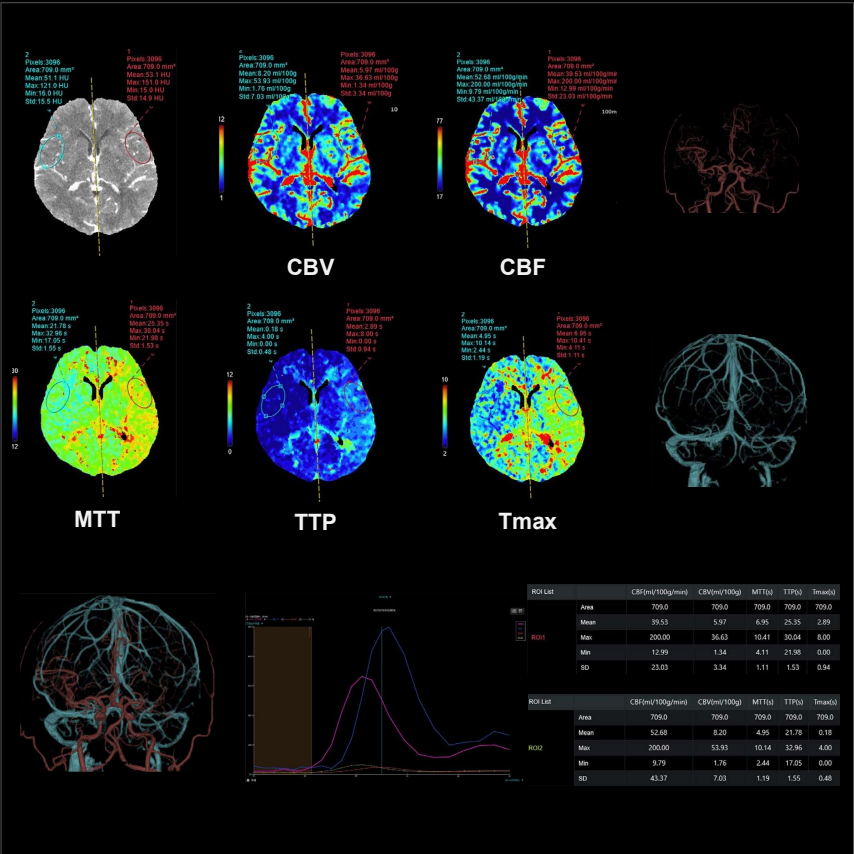

One-stop low-dose whole brain perfusion imaging

Intelligent automatic acquisition of multiple perfusion parameters including CBV, CBF, TTP, MTT, Tmax in one click.

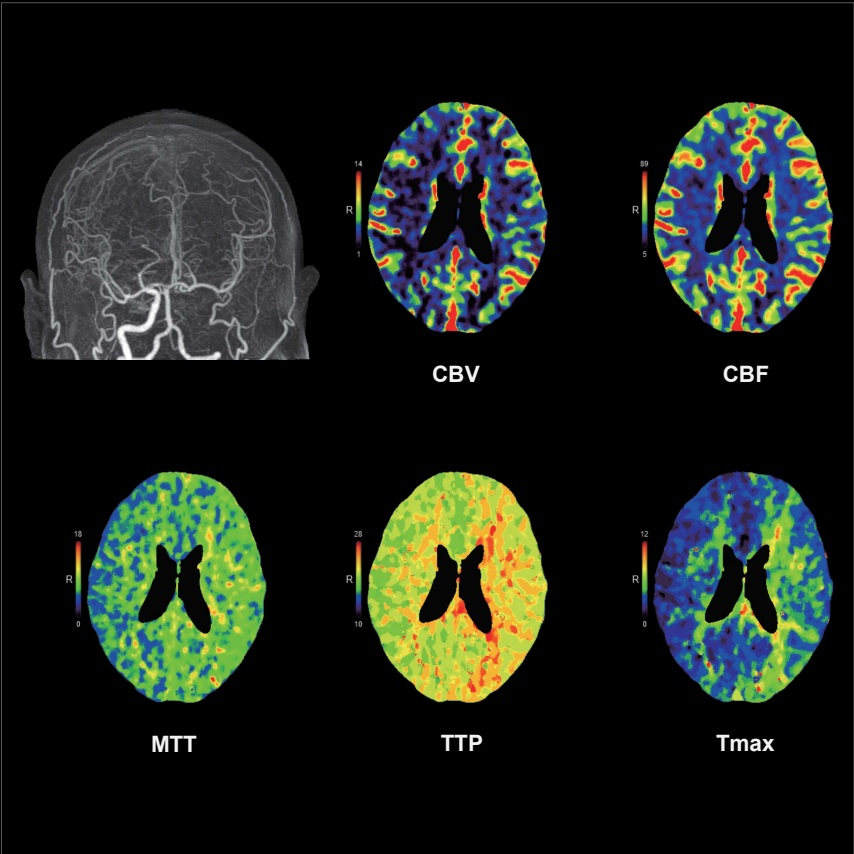

Oncology

Dual energy, Stationary perfusion, Dynamic perfusion and more precise diagnosis.



Dual-Energy Scan

With a 16 cm detector, provide information from anatomy to function and boost post processing for faster and more precise diagnosis.



Stationary Perfusion Scan

The uCT 960+ enables a whole organ perfusion acquisition with up to 160 mm coverage without moving the table. This allows for a perfusion acquisition of the brain, liver, kidneys, pancreas and other organs with uniform contrast along the z-direction.



Dynamic Perfusion Scan


For a scan range over 160 mm, the uCT 960+ enables dynamic perfusion scanning in the helical mode with the table moving back and forth at a variable pitch and supports the a maximum 40 cm scanning range with this scan mode.



Oncology


Dual-Energy Scan

CNR curve can be generated by drawing ROI, and also generate scatterplot, histogram, curves to determined the consistency of the tumor.

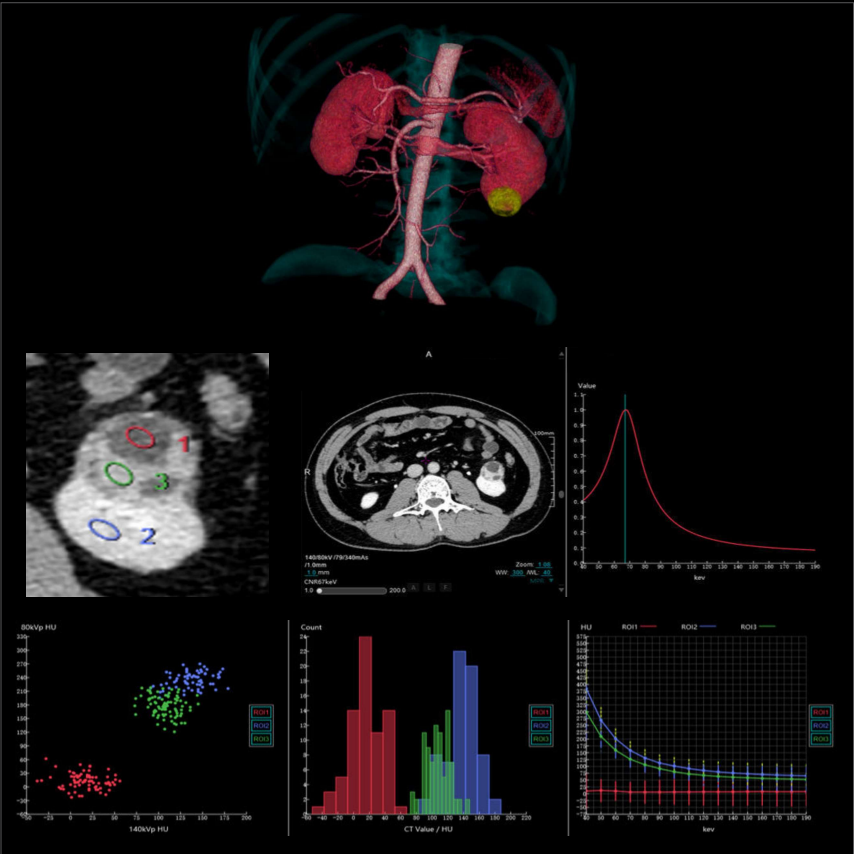

Dual-Energy Scan

The curves clarify the different nature of the bilateral renal lesions. The changes in iodine concentration demonstrated that right renal cancer and left renal cyst.

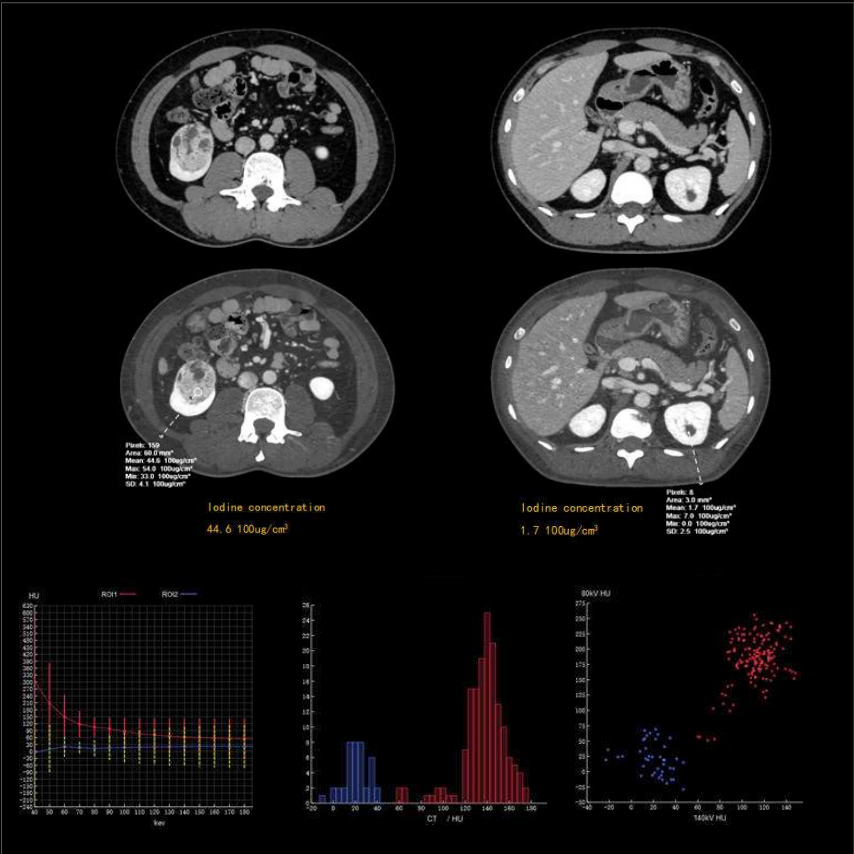

Emergency & Trauma

The uCT 960+ has a capability of fast whole body scans, providing imaging solutions with fast and multi-dimensional information for the clinical decision support in the ED.



Single Rotation Brain Imaging

The entire scan can be completed in a half-second to minimize the motion artifacts caused by involuntary movement of trauma, geriatric or other patients.



One-Stop Acquisition For Triple-rule-out Examination

Using ECG-gated axial scans and non-ECG-gated axial or helical imaging with a fast switching time with only one injection of contrast media. One-stop acquisition for the evaluation of coronary arteries, pulmonary artery and aorta for triple-rule-out examination.




Real Time 3D

Real Time 3D is an innovative solution that will automatically produce MPR or VR preview images in real-time (along with the axial real-time preview images) as the acquisition is taking place.



Intelligent and Precise Analysis


With comprehensive and accurate Intelligent rib fracture analysis, the advanced application platform greatly improves the efficiency of post-processing.

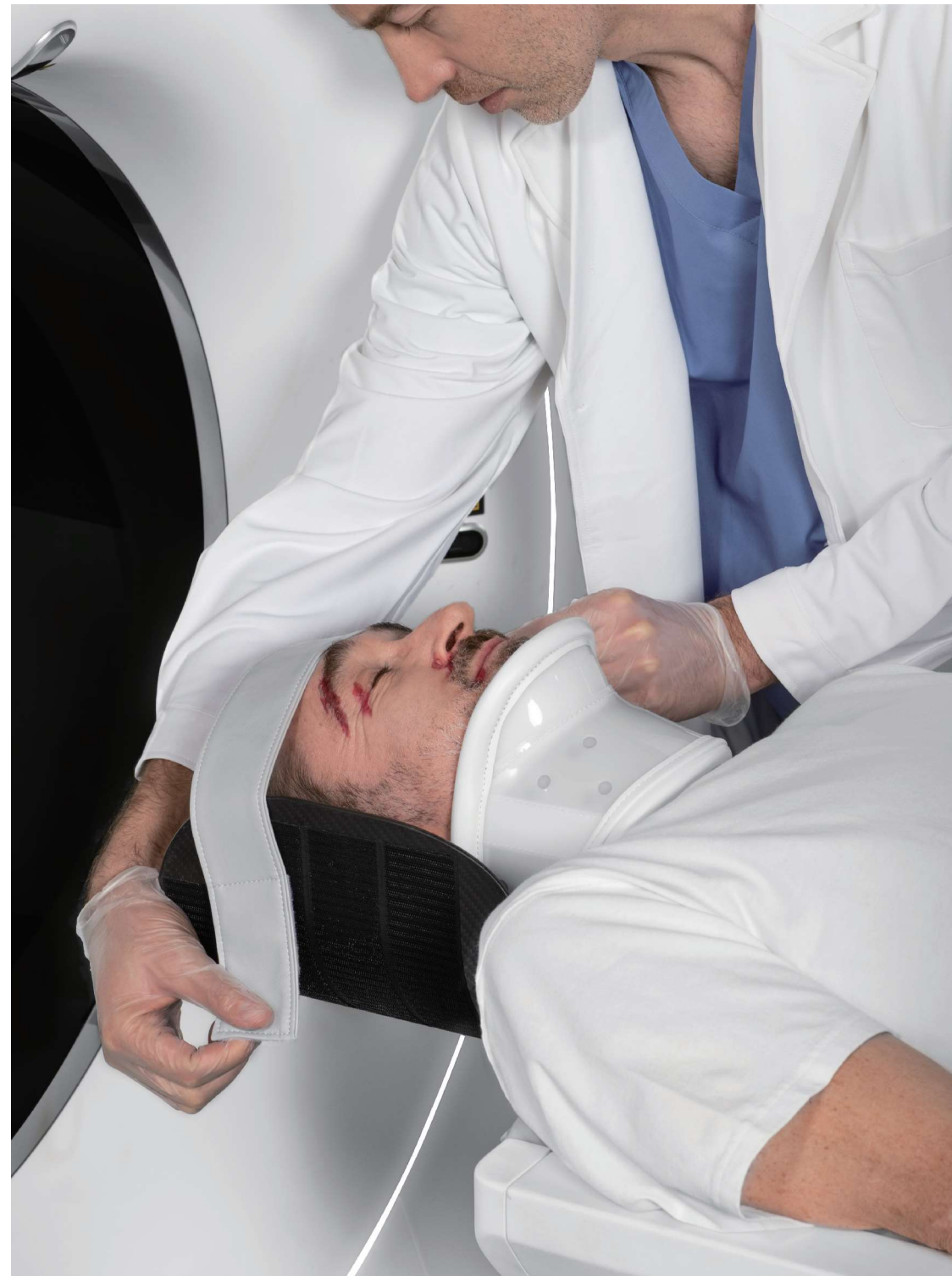

Emergency & Trauma


Brain scan in half a second

16cm whole brain coverage;
Low scattering artifacts at 44.8mGy dose, clearly showing bleeding lesions.

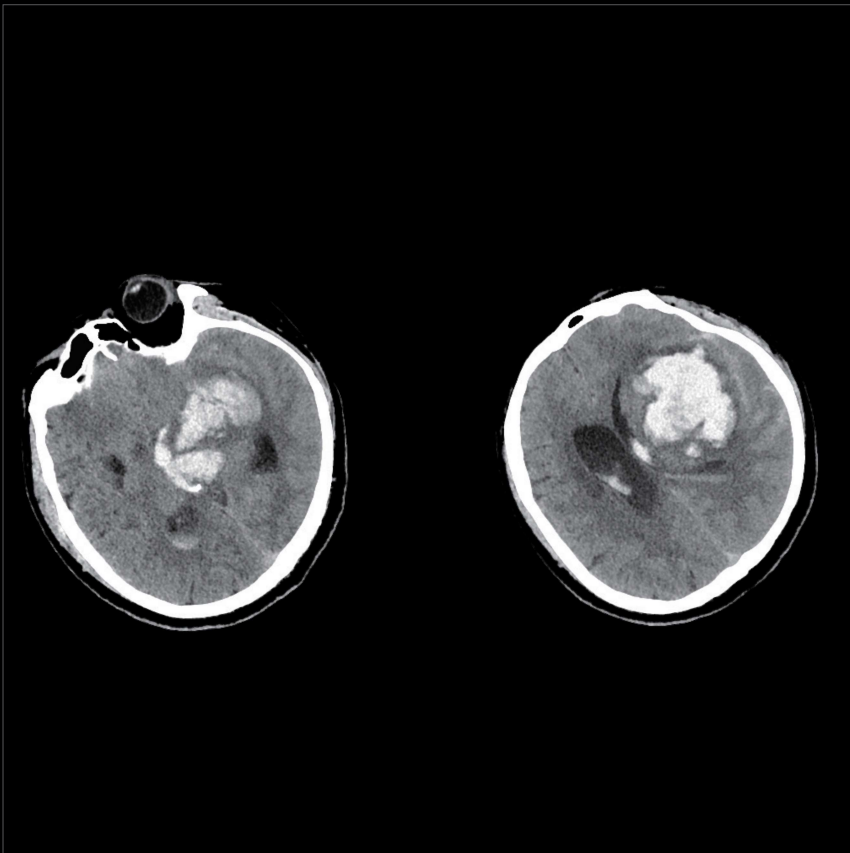

Chest scan in one second

Automatic cutting, opening, straightening and automatic marking of ribs and boost post processing for faster and more precise diagnosis to diagnostic confidence.

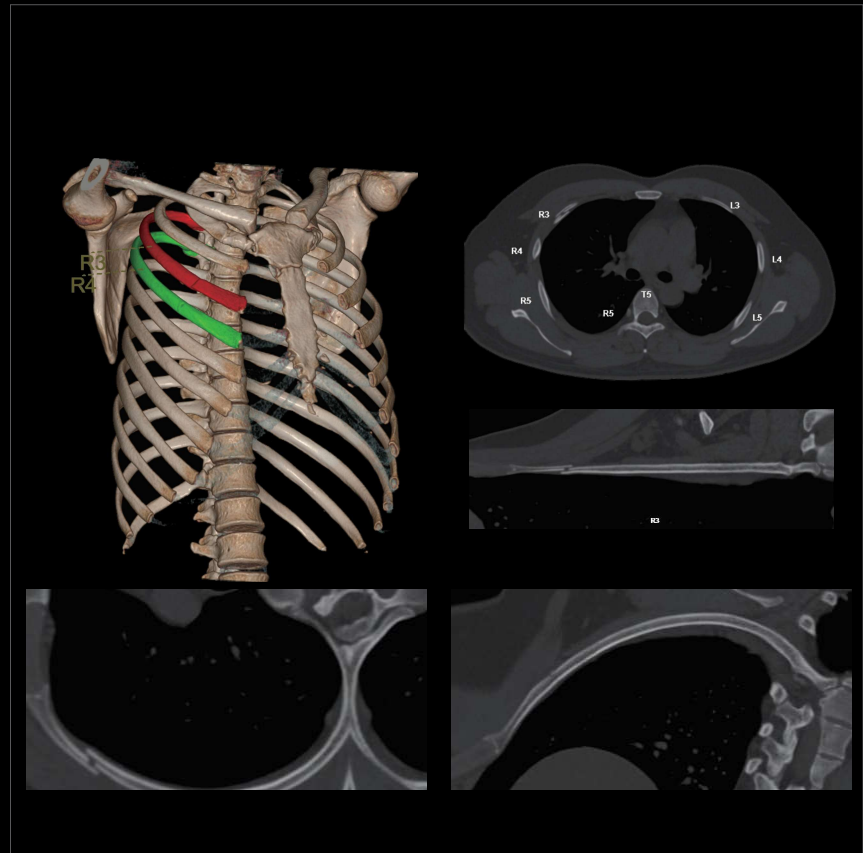

Emergency & Trauma

Abdomen scan in 2 seconds

The 440 mm/s travel speed, this facilitates the scanning of large patients, enables fast scanning for large anatomic ranges.



Large-scale scan in 3 seconds

3 seconds helical scan, quickly obtain high quality images.





First-class experience

The wide 82 cm bore provides a more comfortable examination experience, and flexible operation space.



Low Dose Technology

60 kVp scan capabilities reduce radiation dose while maintaining image quality especially for smaller adults and pediatric patients.




High patient throughput

With the advanced liquid-metal bearing technology the uCT 960+ X-ray tube allows for efficient heat dissipation and reliable performance the X-ray tube of the uCT 960+ has the equivalent anode heat capacity of 30 MHU, and ensures fluent scanning with high patient throughput and consecutive high power examinations.




First-class Experience

60 kV scans allow you to lower the radiation dose while improving image quality and diagnostic confidence.

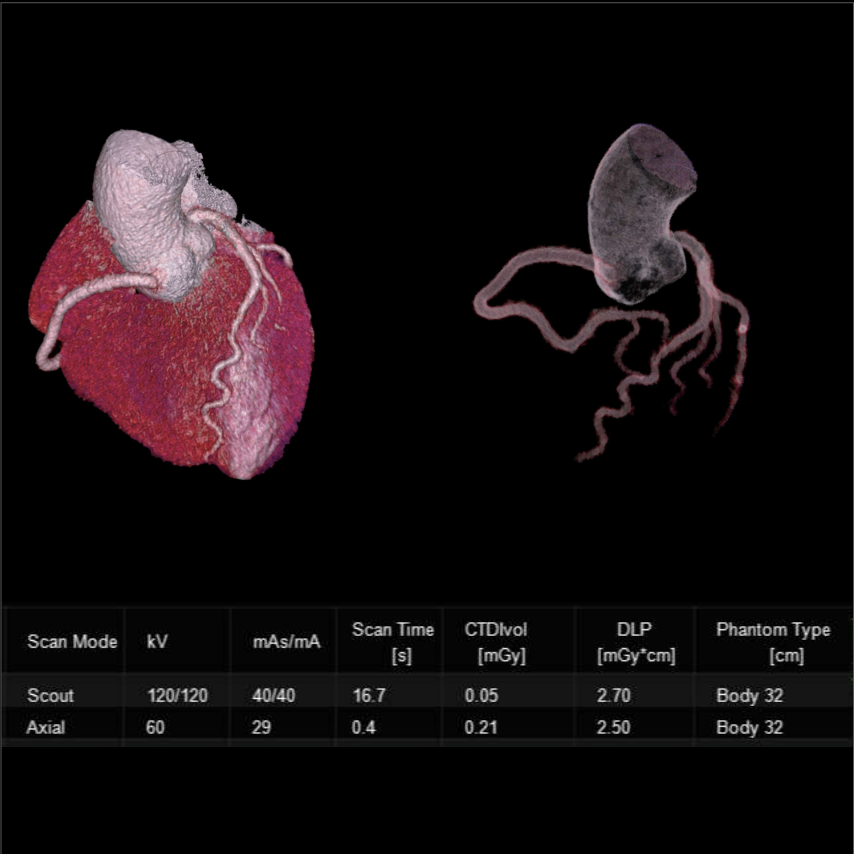


152 kg , BMI 48,82 cm bore size, 300 kg maximum table load capacity allows flexible positioning and access for fat patients.

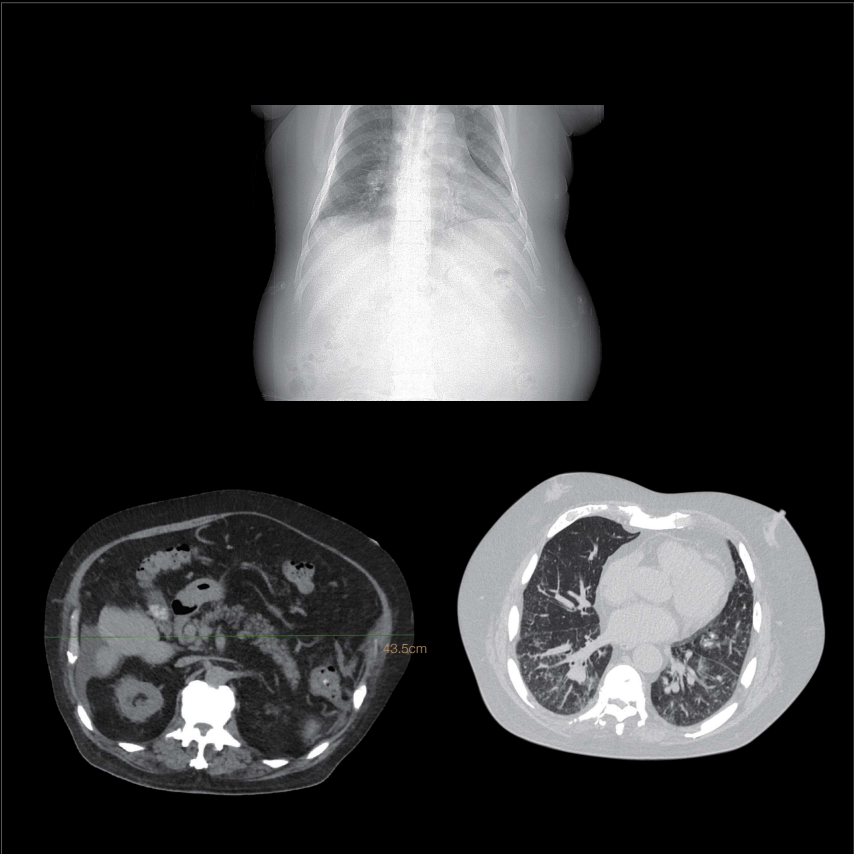

Human-Centered Design

Focusing on user experience, uCT 960+ combines accurate operation with a lightweight and artistic design. We bring aesthetic enjoyment and ease of use to the technology, delivering care, trust and respect through our design.

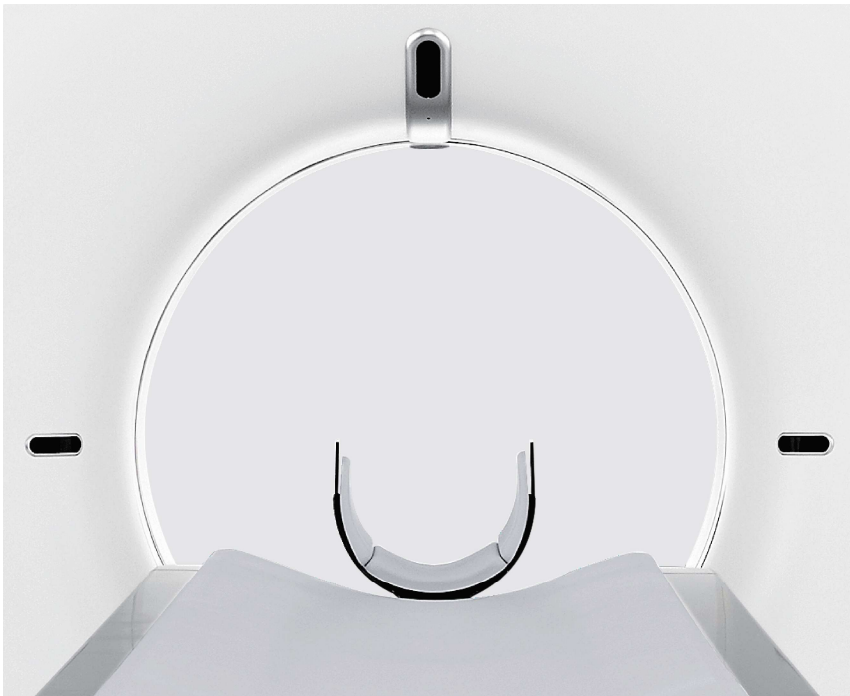

Pleasing Aesthetics

Our design scheme integrates oriental aesthetics with minimalism, presenting a seamless fusion of traditional and modern styling.

User-Friendly Design

The product design delivers comfort, safety, efficiency and ease-of-use. By applying ergonomic principles the uCT 960+ combines innovative design with optimal functionality in order to provide a better possible user experience, optimizing patient comfort during the examination.

Sophisticated Craftsmanship

Driven by the tenets of precision design, we fine-tune every technological detail to embody the spirit of craftsmanship in every product.